

# Bridging Signal and Intelligence: A DSP-Driven EEG Feature Pipeline for Neural Network Seizure Detection

## Presented By

Md. Ohiduzzaman Sajol (2003043)

Durjoy Sarkar Dhrubo (2003026)

## Presented On

November 2, 2025

## Presented To

Dr. Monira Islam

*Associate Professor, Dept. of EEE, KUET*

Dr. Md. Mahabub Hasan

*Professor, Dept. of EEE, KUET*



# Outline

- **Abstract**
- **Introduction**
- **Background**
- **Methodology**
- **Results**
- **Discussion**
- **Conclusion**
- **Future Work**
- **References**

# Abstract

- **DSP-based EEG pipeline for seizure detection**
- **CHB-MIT dataset with notch & bandpass filters**
- **Extracted spectral & statistical EEG features**
- **Visualized raw vs. filtered & frequency spectra**
- **Features ready for ML seizure classification**

# Introduction

- EEG is a **non-invasive** method to measure brain activity.
- Automated EEG-based **seizure detection** is a key biomedical DSP task.
- EEG's **non-stationary** nature demands **strong preprocessing & feature extraction**.
- DSP provides tools for **time–frequency** domain analysis (e.g., **convolution, DFT, filtering**).
- EEG data (**CHB-MIT**) is preprocessed, features extracted, and prepared for ML seizure classification.

# Background: EEG and Seizure Activity

- EEG records brain's electrical signals via scalp electrodes.
- Typical EEG range: **0.5–70 Hz**, divided into standard bands:
  - $\Delta$  (**0.5–4 Hz**) – Deep sleep
  - $\theta$  (**4–8 Hz**) – Drowsiness/meditation
  - $\alpha$  (**8–13 Hz**) – Relaxed wakefulness
  - $\beta$  (**13–30 Hz**) – Active focus
  - $\gamma$  (**30–70 Hz**) – Cognitive processing
- **Seizures:** sudden, synchronized bursts of activity
  - visible as amplitude spikes & spectral power changes

# Methodology: Data Acquisition

- Dataset: **CHB-MIT Scalp EEG Database** (pediatric epilepsy patients)
- Format: **EDF files** with multi-channel EEG recordings
- Loaded via *mne.io.read\_raw\_edf()*
- Includes both **seizure** and **non-seizure** intervals

# Methodology: Preprocessing Pipeline

- **Notch Filter (60 Hz):** remove power-line noise (IIR notch)
- **Bandpass FIR Filter (0.5–70 Hz):** isolate standard EEG range
- **Segmentation:** 10-s overlapping epochs for temporal analysis
- **Normalization:** z-score (zero mean, unit variance) for consistency

# Methodology: DSP Operations

- **Correlation:** inter-channel Pearson correlation matrices
- **Convolution:** moving-average smoothing to reduce fluctuations
- **Spectral Analysis:** FFT + Welch PSD for frequency and band-power estimation



# Methodology: Feature Extraction

- **Band Powers:** mean, std, median in  $\delta$ ,  $\theta$ ,  $\alpha$ ,  $\beta$ ,  $\gamma$  bands
- **Statistical:** energy, variance, skewness, kurtosis
- **Spectral Entropy:** Shannon entropy of normalized PSD
- **Correlation Metrics:** mean, std, max of inter-channel correlations
- **Band Ratios:** e.g.,  $\alpha / \beta$  for cognitive state insight


# Methodology: Implementation

- Tools: **Python 3.10**, MNE, SciPy, NumPy, pandas, Matplotlib
- Interactive workflow via **Jupyter Notebook**
- Features aggregated into *chb\_features.csv*
- Repository: [github.com/projectohid/4142-project](https://github.com/projectohid/4142-project)

## Results: Spectrogram Analysis (Normal vs. Seizure)


- **Normal EEG (Fig. 1)** : Broad spectral power (10–30 Hz,  $\alpha$  &  $\beta$  bands)
- **Seizure EEG (Fig. 2)** : Dominant low-frequency power (<10 Hz,  $\delta$  &  $\theta$ )
- Confirms effective filtering and **distinct physiological vs. pathological patterns**

# Results: Spectrogram Analysis (Normal vs. Seizure)

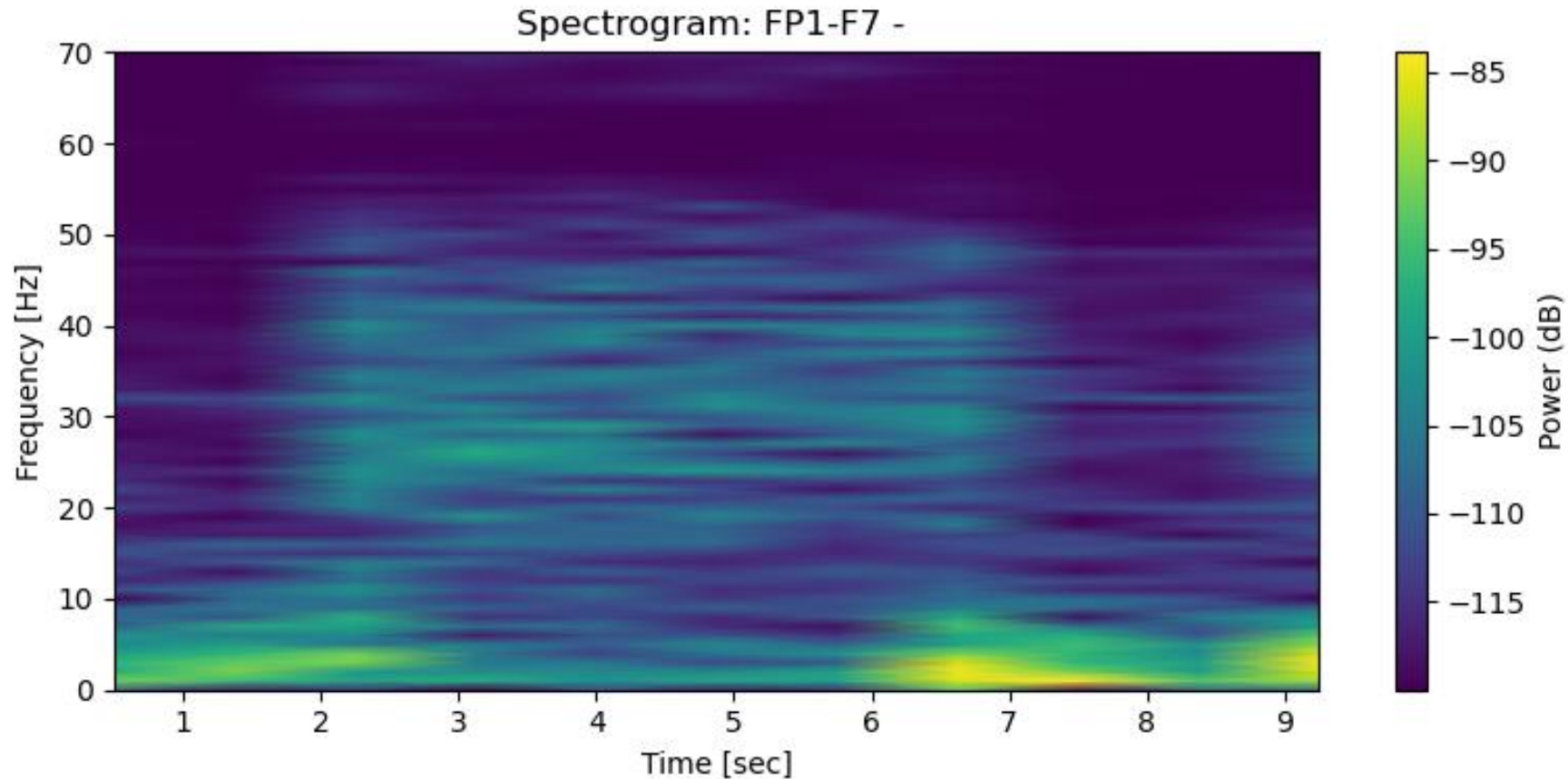

**Fig. 1: Spectrogram for *Chb01\_01.edf* (non-seizure).**

# Results: Spectrogram Analysis (Normal vs. Seizure)

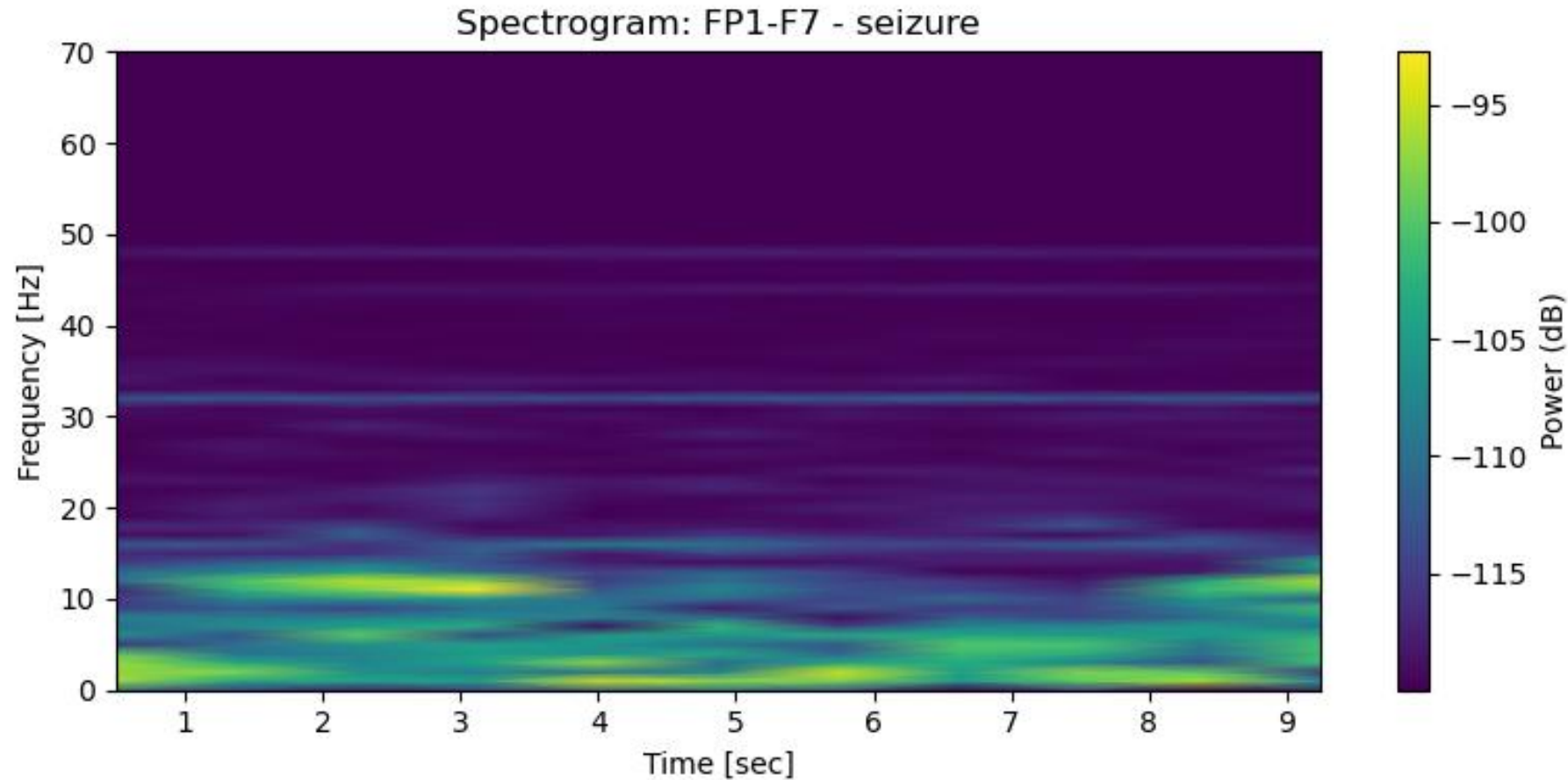


**Fig. 2: Spectrogram for *Chb01\_03.edf* (seizure).**

## Results: Raw Vs. Filtered Signal

- **Time-domain (Fig. 3):** Filtered signal smoother; reduced drift & 60 Hz noise
- **Frequency-domain (Fig. 3):** 60 Hz peak suppressed; spectrum confined to 0.5–70 Hz
- Verifies **notch + FIR bandpass filters** effectiveness

# Results: Raw Vs. Filtered Signal

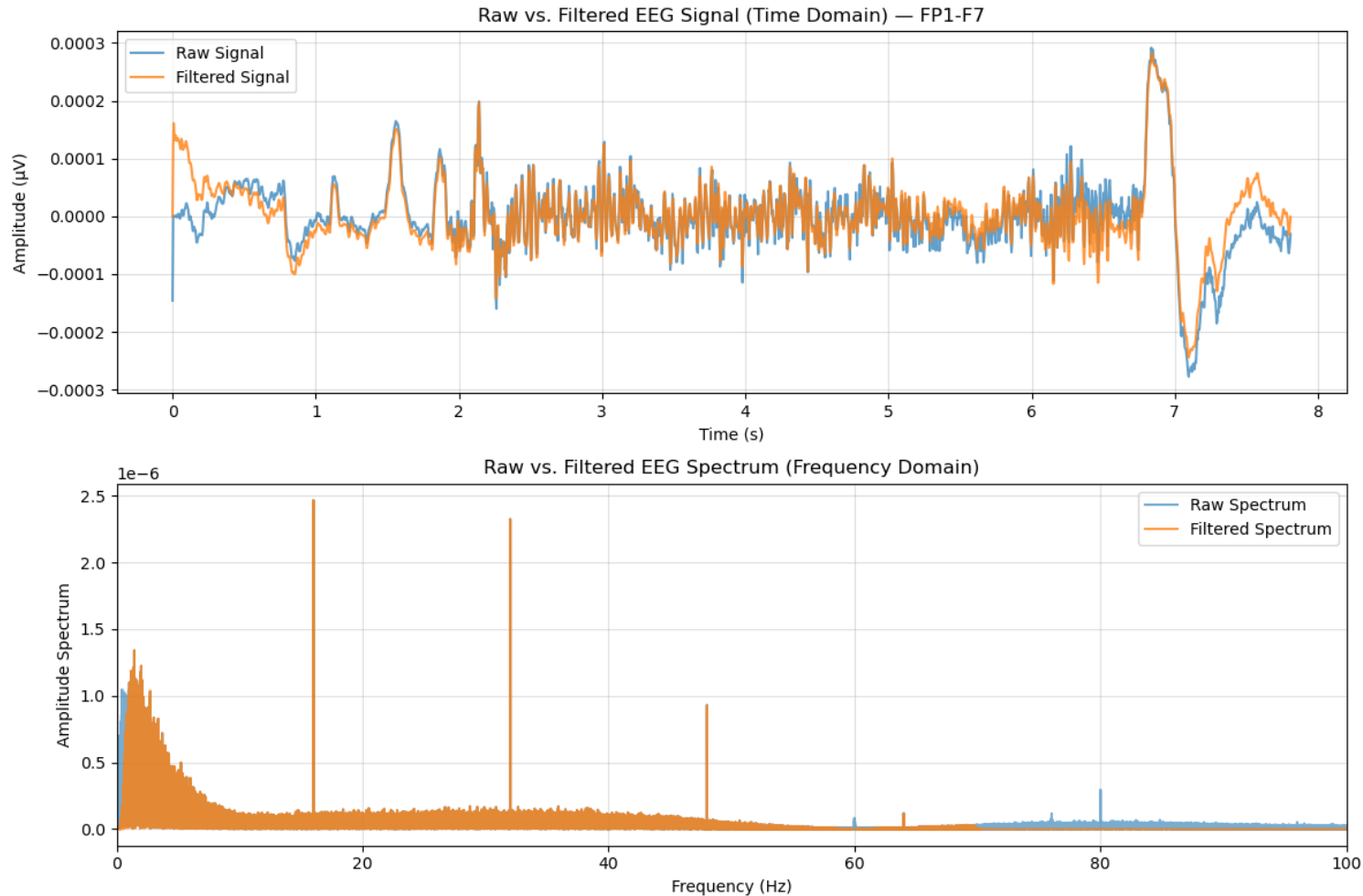


**Fig. 3: Raw and filtered Signals and their Spectrums for *Chb01\_01.edf*.**

# Results: EEG Band Power Distribution

- Power hierarchy:

$$\Delta > \theta > \beta > \gamma > \alpha$$

- Low-frequency dominance  
→ seizure or drowsy state indicator
- Confirms **spectral accuracy**  
of feature extraction

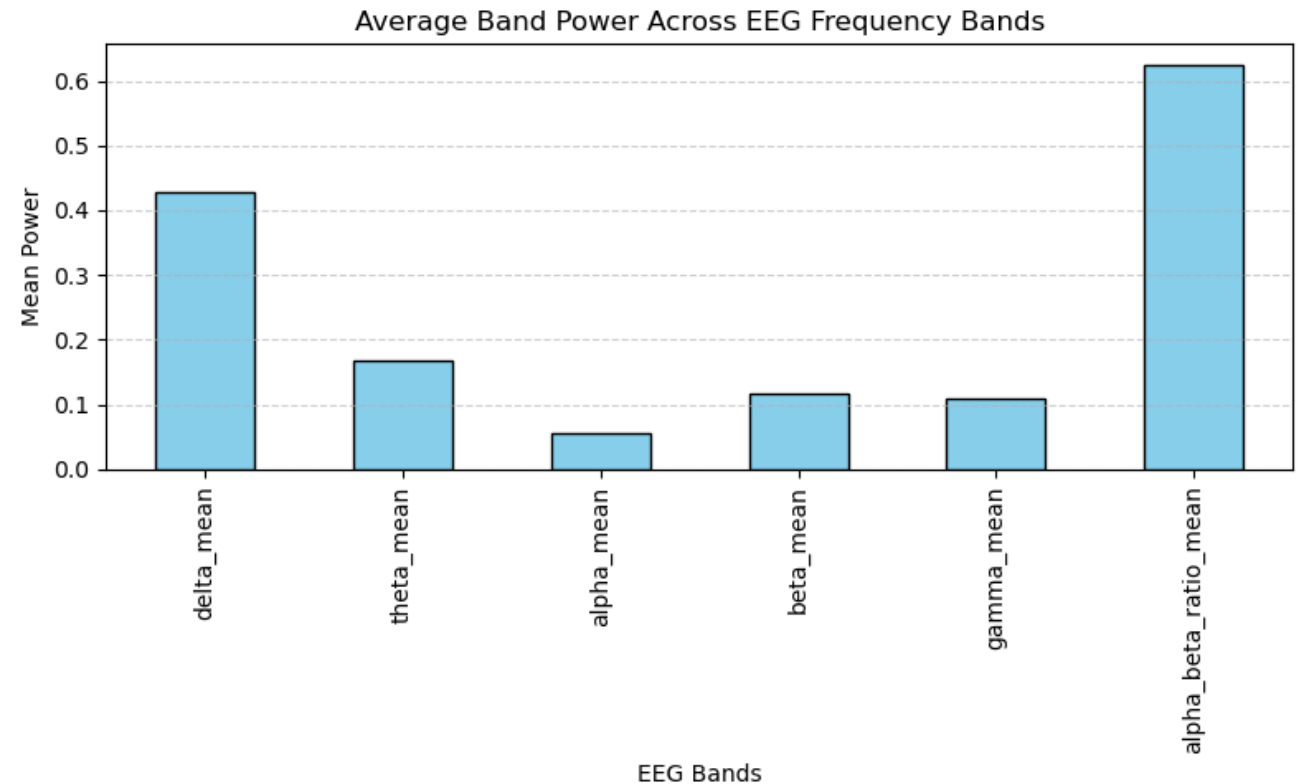


**Fig. 4: Average Band Power Across EEG Frequency Bands for *Chb01\_01.edf*.**



# Results: Time-Domain Statistical Features


- **High kurtosis:** sharp, transient seizure spikes
- **Low skewness:** symmetric yet with frequent outliers
- Confirms **impulsive, non-Gaussian EEG** patterns in seizures

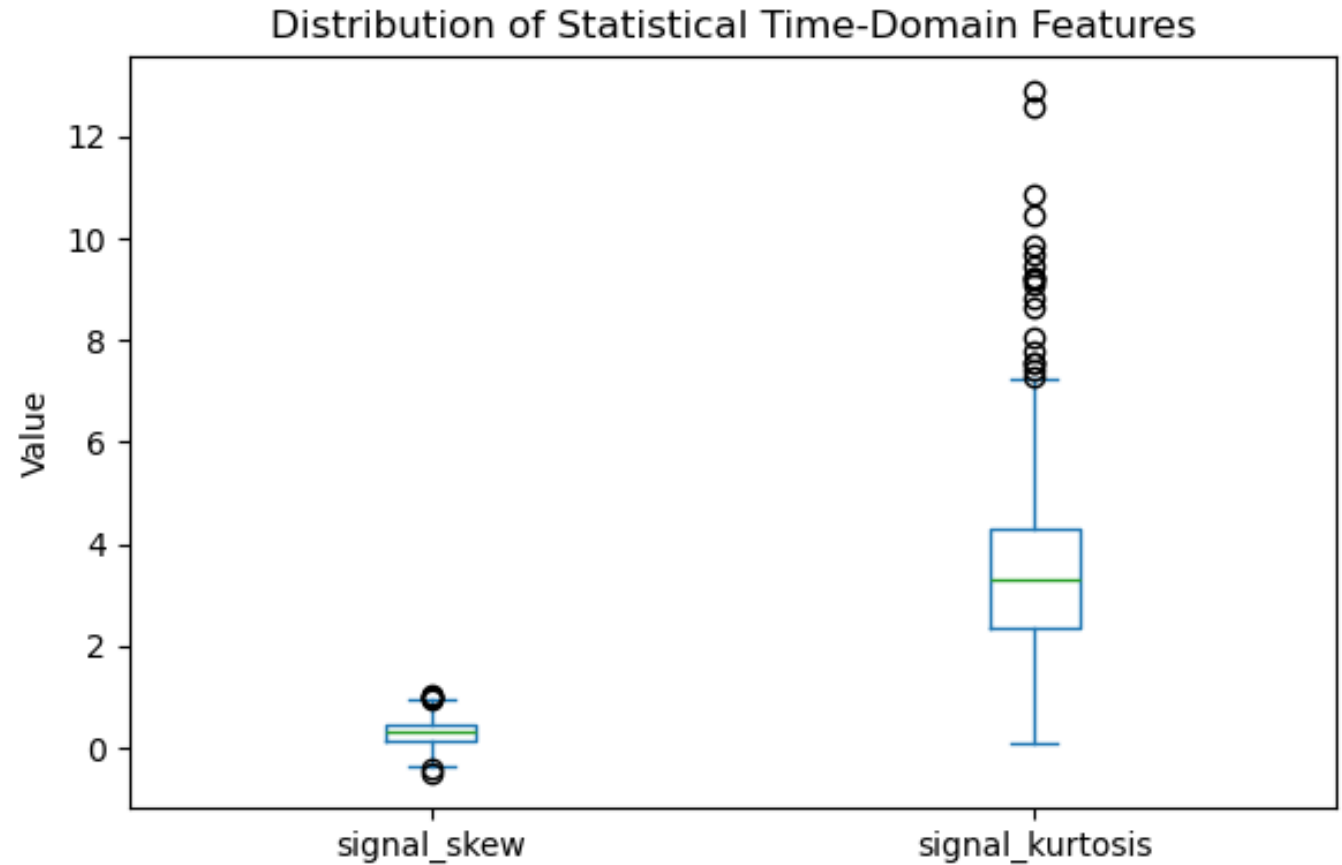

**Fig. 5: Distribution of Statistical Time-Domain Features for *Chb01\_01.edf*.**

# Results: Spectral Entropy vs. Delta Power


- **Negative correlation** → higher  $\delta$  power = lower entropy
- Indicates **reduced complexity** during seizure activity
- Validates **entropy as nonlinear brain state marker**

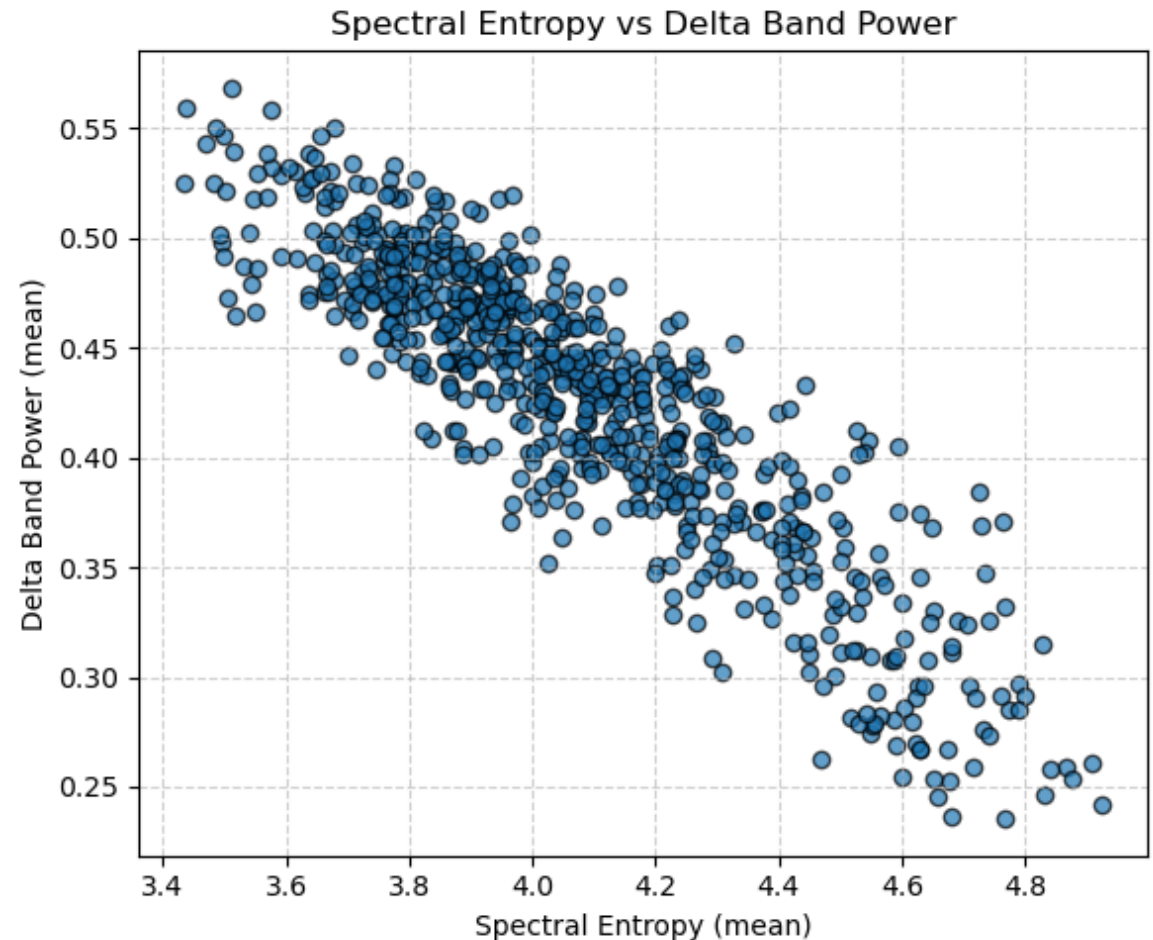

**Fig. 6: Spectral Entropy vs. Delta Band Power for *Chb01\_01.edf*.**

# Results: Inter-Channel Correlation

- **Max correlation  $\approx 1.0$ :**  
strong cross-channel  
synchronization
- **Consistent** mean correlation  
= stable phase coupling
- Demonstrates **spatial  
coherence** in seizure EEG

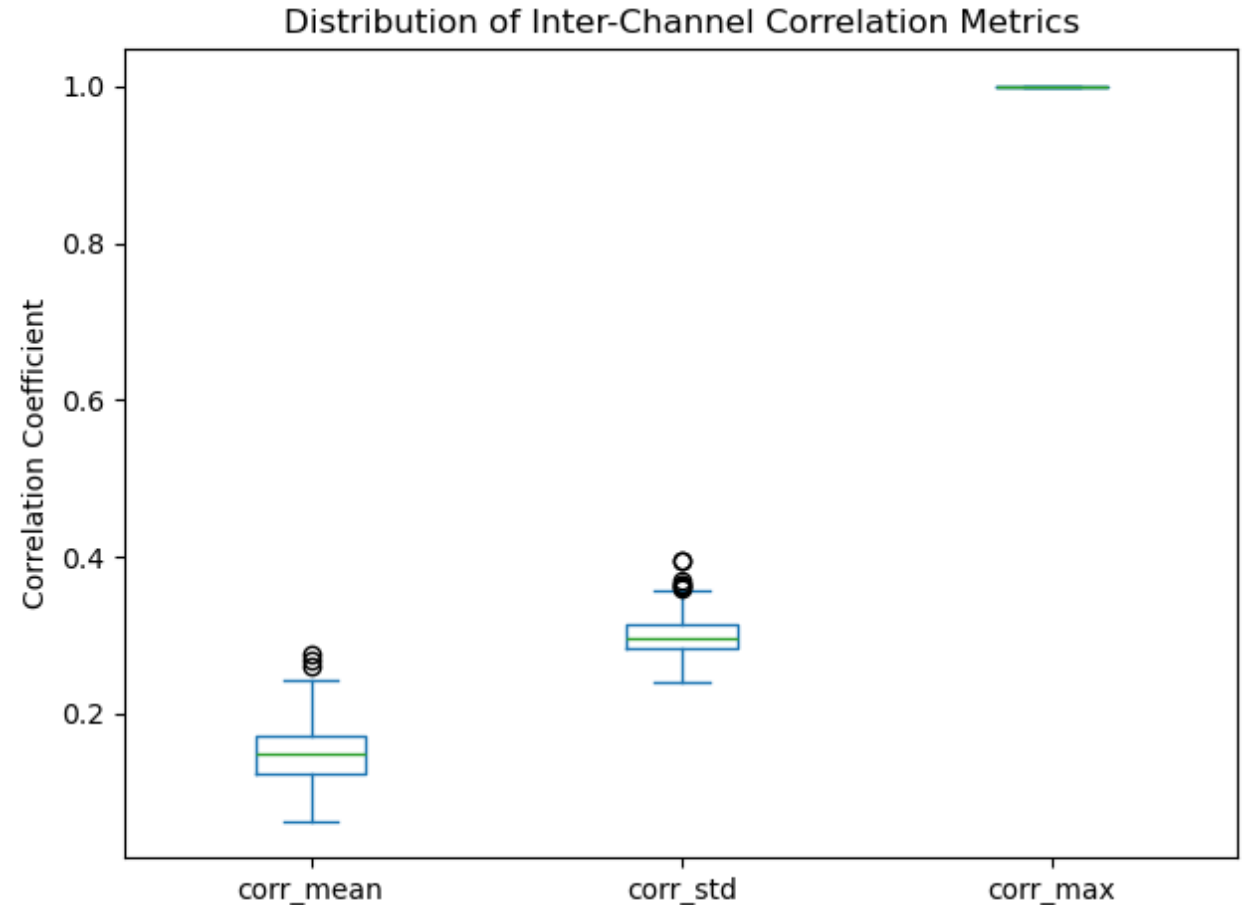

**Fig. 7: Distribution of Inter-Channel Correlation Metrics for *Chb01\_01.edf*.**

# Results: Feature Matrix & ML Readiness

**Table I: Part of the extracted Feature Matrix for *Chb01\_01.edf*.**

	delta_mean	delta_std	delta_max	delta_median	theta_mean	theta_std	theta_max	theta_median	alpha_mean	alpha_std
0	0.439917	0.122878	0.680715	0.425868	0.134942	0.048450	0.229849	0.132495	0.046567	0.029513
1	0.491996	0.147472	0.753315	0.489125	0.157016	0.071678	0.355560	0.139707	0.052553	0.036122
2	0.375810	0.082527	0.485999	0.389008	0.155521	0.067377	0.320863	0.145871	0.055865	0.031640
3	0.328070	0.089442	0.475458	0.338504	0.164152	0.065626	0.293830	0.164945	0.056487	0.021148
4	0.467772	0.085527	0.657703	0.458547	0.173338	0.052196	0.284436	0.168564	0.054683	0.029469

5 rows × 34 columns

- Each 10 s epoch → **30+ DSP-based features**
- Seizure/non-seizure epochs clearly separable in feature space
- Confirms **pipeline suitability for neural network models**

# Results: Summary

Domain	Key Observation	DSP Technique
Spectrograms	Low-frequency dominance	STFT / Spectrogram
Filtering	Removed 60 Hz & artifacts	Notch + FIR
Band Power	High $\delta$ , $\theta$ ; low $\alpha$ , $\beta$ , $\gamma$	FFT + PSD
Statistics	High kurtosis, low skewness	Time-domain
Entropy	–ve correlation with $\delta$ power	Entropy calc.
Correlation	Strong inter-channel sync	Pearson corr.
Feature Set	Clear seizure separation	DSP feature extraction

# Discussion

- **Spectral Behavior**

- Normal EEG: broad power (10–30 Hz,  $\alpha$ – $\beta$  bands)
- Seizure EEG: strong low-freq. dominance (<10 Hz,  $\delta$ – $\theta$ )
- Confirms **shift from desynchronized** → **synchronized** neuronal activity
- Filtering (notch + bandpass) effectively removed 60 Hz noise & drift while preserving core EEG content

- **Band Power Trends**

- Hierarchy:  $\Delta > \theta > \beta > \gamma > \alpha$
- High  $\alpha/\beta$  ratio & reduced high-freq. power → decreased cortical processing during seizures

# Discussion

- **Time-Domain Statistics**
  - **High kurtosis:** sharp epileptic transients
  - **Low skewness:** symmetric but heavy-tailed amplitude distribution
  - Indicates **impulsive, non-Gaussian** seizure morphology
- **Entropy & Complexity**
  - **Entropy**  $\downarrow$  as  **$\Delta$ -power**  $\uparrow \rightarrow$  reduced signal complexity during seizures
  - Captures transition to **predictable, synchronized** brain dynamics
  - Entropy = nonlinear complement to power-based metrics

# Discussion

- **Inter-Channel Correlation**
  - **High cross-channel correlation ( $\approx 1.0$ )** → strong spatial coherence
  - Reflects widespread cortical synchronization in seizure states
- **Integrated DSP Takeways**
  - **Filtering:** ensures clean, physiologically valid signals
  - **FFT / PSD:** reveal spectral composition
  - **Statistics + Entropy:** quantify morphology & complexity
  - **Correlation:** captures spatial coupling → Together, **classical DSP techniques uncover physiological + pathological EEG patterns**, producing interpretable features for intelligent seizure detection.



# Conclusion

- DSP methods effectively applied for EEG seizure analysis.
- Filtering + spectral + statistical tools made EEG **interpretable**.
- **Key Points:**
  - **Preprocessing:** Notch & FIR filters removed noise, kept key signals.
  - **Spectral:** Clear normal–seizure separation (low-freq dominance).
  - **Stats:** High kurtosis, low entropy → synchronized seizure waves.
  - **Correlation:** Strong inter-channel coupling during seizures.
  - **Features:** Compact, ML-ready dataset for classification.
- **Takeaway:**
  - DSP bridges theory & application, enabling explainable EEG systems.

# Future Work

- **Wavelet Filters:** Handle dynamic, non-stationary EEG noise.
- **ICA Cleaning:** Remove eye/muscle artifacts.
- **Deep Learning:** Train models on DSP-derived features.
- **Automation:** Streamline preprocessing & feature selection.

# References

- [1] L. Tan and J. Jiang, *Digital Signal Processing: Fundamentals and Applications*, 3rd ed. Elsevier Science, 2018.
- [2] M. Teplan, “Fundamentals of EEG measurement,” *Measurement Science Review*, vol. 2, no. 2, pp. 1–11, 2002.
- [3] H. Adeli, Z. Zhou, and N. Dadmehr, “Analysis of EEG records in an epileptic patient using wavelet transform,” *Journal of Neuroscience Methods*, vol. 123, no. 1, pp. 69–87, 2003.
- [4] U. R. Acharya, S. V. Sree, G. Swapna, R. J. Martis, and J. S. Suri, “Automated EEG analysis of epilepsy: A review,” *Knowledge-Based Systems*, vol. 45, pp. 147–165, 2013.
- [5] A. Aarabi, R. Fazel-Rezai, and Y. Aghakhani, “EEG seizure prediction: Measures and challenges,” *Computers in Biology and Medicine*, vol. 42, no. 4, pp. 376–390, 2012.

# Thank You

## Any Questions?

